

STANDUP Series

AV 11 1110 STAND UP 1110 gas spring height adjustment, manual lift

This device enables the user to assume and maintain an upright position, and to move indoors assisted by a carer.

The lift is achieved manually, which means that the user can grip the side handles and the tray edges.

It features large handles as well as a wide tray to hold objects¹.

If the user is unable to bring themselves to a vertical standing position, they must be helped by one or more carers depending on their body weight. A series of adjustments facilitate accessibility² and allow the aid to be adapted to the user's morphological characteristics³.

In particular, the main body can adjust to the user's height by means of a gas spring assisted system⁴. For further information please refer to the presentation on general characteristics, adjustments and technical data.

Weight: 50 kg

Maximum load: 140 kg

User height: min.140 cm - max 200 cm



1) - Wide tray with shaped lateral grips and top recesses to hold objects



2) - The sliding footplates move towards the user when transferring from the wheelchair



3) - Telescopic tray adjustment system (height adjustment).




4) - Gas spring control to adjust the height of the support structure.



POSSIBILITY OF UPGRADING THE MODELS FROM THE STAND UP SERIES

If the user's residual capabilities decrease over time, it is possible to upgrade the model AV 11 1110 to one of the models with servo-assisted lift. The upgrade is easy, but it must be carried out by a skilled technician. The standing frame is already arranged to fit any necessary subsequent additions that can be of great help to the user.

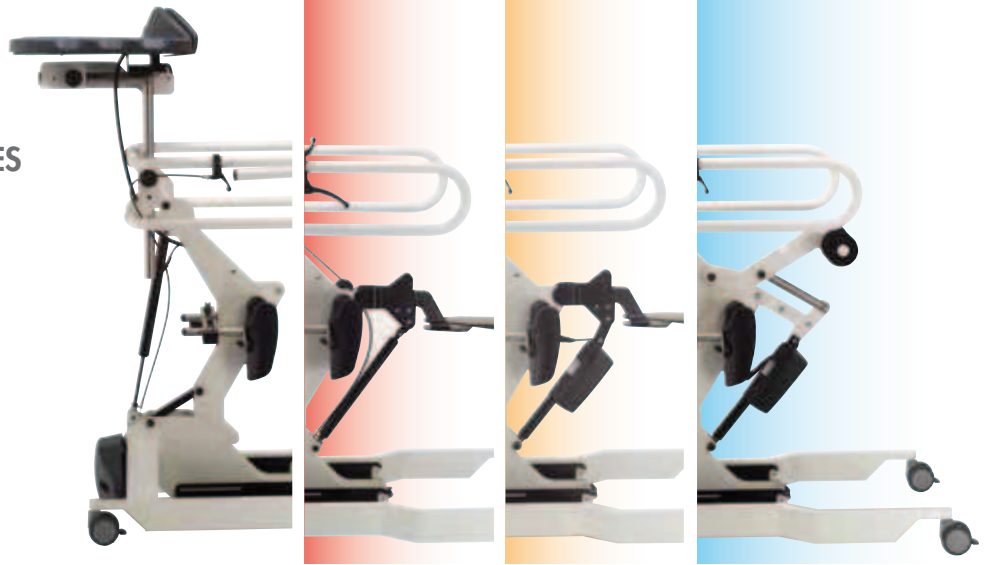



- lift by gas strut and seat
- electric lift with seat
- electric lift with harness

ILLUSTRATION OF THE RESTRAINT SYSTEM

A safe restraint of the user when in upright position is ensured by a lumbosacral support that is closed using an accessible handle and locking system.



TECHNICAL DATA

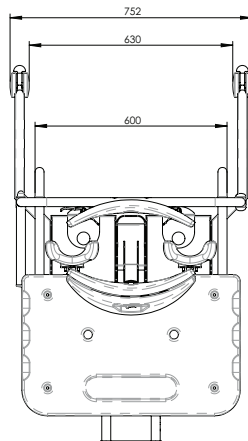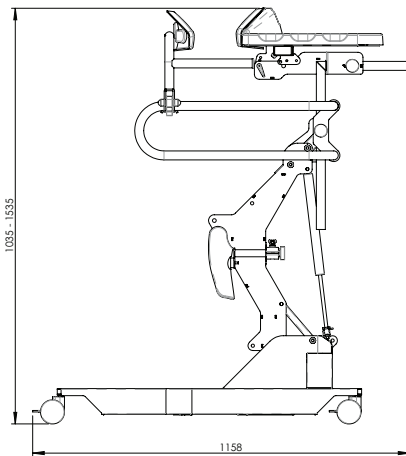

DESCRIPTION	STAND UP SERIES
Lift	manual
Structure height adjustment	gas
Parking brake	on all 4 castors
Frame height from ground	31 mm


ACCESSORIES



AC0014
LEG STRAP
to stabilize the legs



AC0044
BACK STRAP M
to restrain the trunk



AC0047
FOOT STRAP
to stabilize the feet position