

LEVABO[®]

medical cushions

- Limited budget
- Scarce resources

Try Levabo's product range of inflatable single-patient devices

- Hygienic
- Comfortable
- Easy
- Economic



Distributed in Ireland by




Pressure relief of heels.


Prevention and treatment, category 1-4



Name		Size in cm.	Item no.	Application	About the device	Pressure map
Heel up - Short		18 x 32 x 25	70001		Heel protector for mobile and seated users, typically wheelchair. Can be used for both prevention and treatment of pressure sores.	
Heel up - Medium		18 x 39 x 25	70007		Heel protector useful for both prevention and treatment of pressure sores. Typically for users with short legs. Most useful for bed-bound patients. It supports the fore-foot and minimizes risk of foot drop.	
Heel up - Long		18 x 48 x 25	70002		Heel protector useful for both prevention and mitigation of pressure sores. Most useful for bed-bound patients. It supports the front foot and minimizes the risk of foot drop.	
Heel up - Max Medium		18 x 39 x 30	70035		Heel protector useful for both prevention and mitigation of pressure sores. Most useful for bed-bound patients and for users with a long foot. It has an extra air chamber under the foot, which provides optimal support of the forefoot.	
Heel up - Max		18 x 50 x 30	70018		Heel protector useful for both prevention and mitigation of pressure sores. Most useful for bed-bound patients and for users with a long foot. It has an extra air chamber under the foot, which provides optimal support of the forefoot.	
Heel up - AIR Medium		18 x 39 x 25	70036		Heel protector useful for both prevention and mitigation of pressure sores. Most useful for mobile patients. It has an air chamber under the heel, which protects the heel when touching the base.	
Hevel up - CUT Medium		18 x 39 x 25	70037		Heel protector useful for both prevention and treatment of pressure sores. Most useful for diabetic patients. The sail under foot is without air allowing cutout for diabetic ulcers on the feet and toes.	

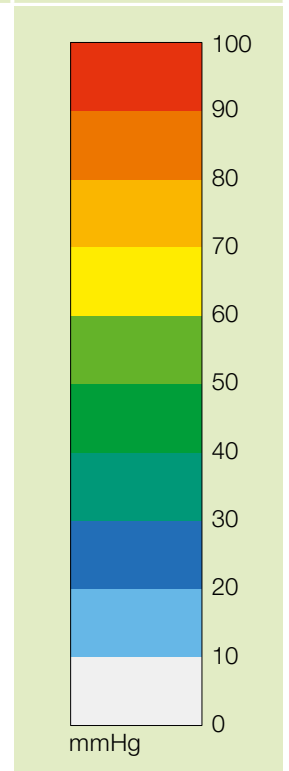

Heel Up® has some important features that provide a number of advantages

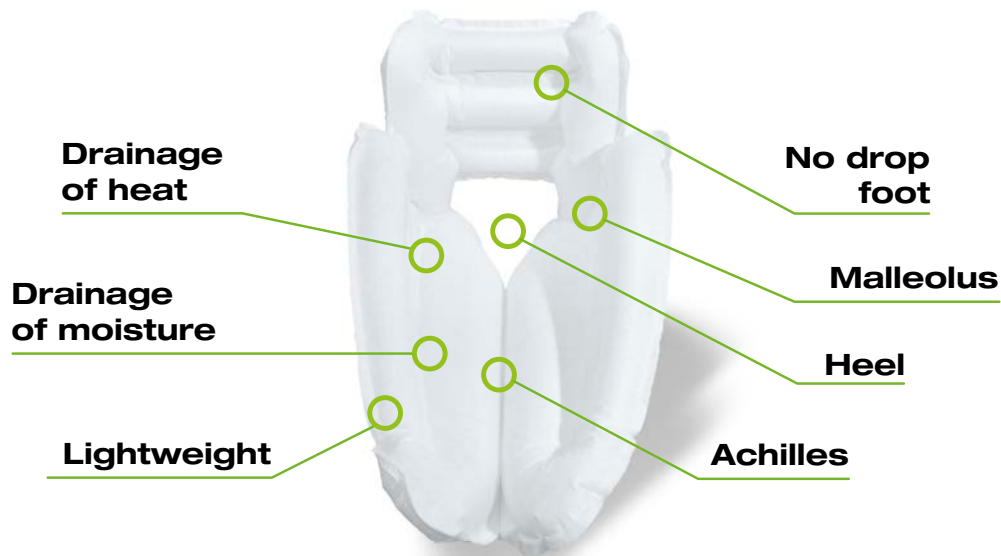

The outer surface provides superior drainage of moisture and heat from the skin, it also ensures the skin doesn't stick to the cushion. It can be wiped with a damp cloth or alcohol.



There is room for the Achilles tendon between the 2 bearing cells. The cells are bevelled so that the lower leg and knee is not unnecessarily burdened.



It is open to the heel and there is room for the malleoli to ensure full relief.



The sail and air under the foot supports the forefoot and minimizes the risk of foot drop.

From zero to pressure relief in less than 60 seconds....



Unfold Heel Up® and insert the pump into the self-closing blue valve until it is tight.



Pump up...



Inflate Heel Up® until the enveloping paper strap - Stop Up® opens - and remove the pump. Heel Up® is now ready to use.



Place Heel Up® on the patient's foot. Close the ribbon across the opening.



Check that the heel is raised about 1½ finger width from the ground.



Done...



Cushions for pressure relief in the seated position





Name	Size in cm.	Item no.	About the device	Application
All up - Seat	42 x 42 x 5	70008	Seat cushion, extra strong. Pressure distribution for the mitigation of pressure sores. Standard wheelchair size. Max 150 kg.	
All up - Seat Large	52 x 52 x 5	70019	Pressure distribution for the prevention and mitigation of pressure sores. Extra long and wide. For patients weighing a maximum of 150kg.	
All up - Seat Max (bariatric)	52 x 62 x 5	70020	Seat cushion, extra strong. Bariatric cushion for user weights up to 250 kg. Pressure distribution for the prevention and mitigation of pressure sores. Extra long and wide.	
All up - Seat Twin	48 x 52 x 5	70021	Seat cushion, extra strong. Suitable for hemiplegic patients with unstable posture. The cushion is split in two to adjust for individual air pressure. For the prevention and relief of pressure sores. Max 150 kg.	
All up - Seat Tailbone	42 x 42 x 5	70022	Special sitting cushion, extra strong. To apply for full relief of the rectum, typically for radiation damage, a sensitive tailbone or great discomfort caused by hemorrhoids. Max 150 kg.	
All up - Back & Seat	27 x 45 x 5	70009	Sitting cushion and back cushion. With advantage, usefull for both areas. Good back support, easy cushion for correct seating. Max 150 kg.	
All up - Donut	45 x 45 x 8	70017	Seat cushion, extra strong. Annular (ring shaped) with no hole in the middle, but with optimum pressure for the relief and prevention of pressure sores. Max 150 kg.	
All up - Bath	50 x 50 x 8	70013	Special seat cushion, extra strong. For bath and toileting chairs to offer a better distribution of pressure for prolonged use. Can be fastened to the chair with tape. After bathing, the cushion can be dried and reused. Max 150 kg.	

Cushions for special positioning

All up - Belly	65 x 75 x 15	70023	Special positioning cushions with particular relevance for radiology departments to elevate the chest during examinations.	
All up - Kidney	50 x 20 x 20	70014	Special positioning cushion with particular relevance for radiology departments for supporting the kidney/pelvic area.	




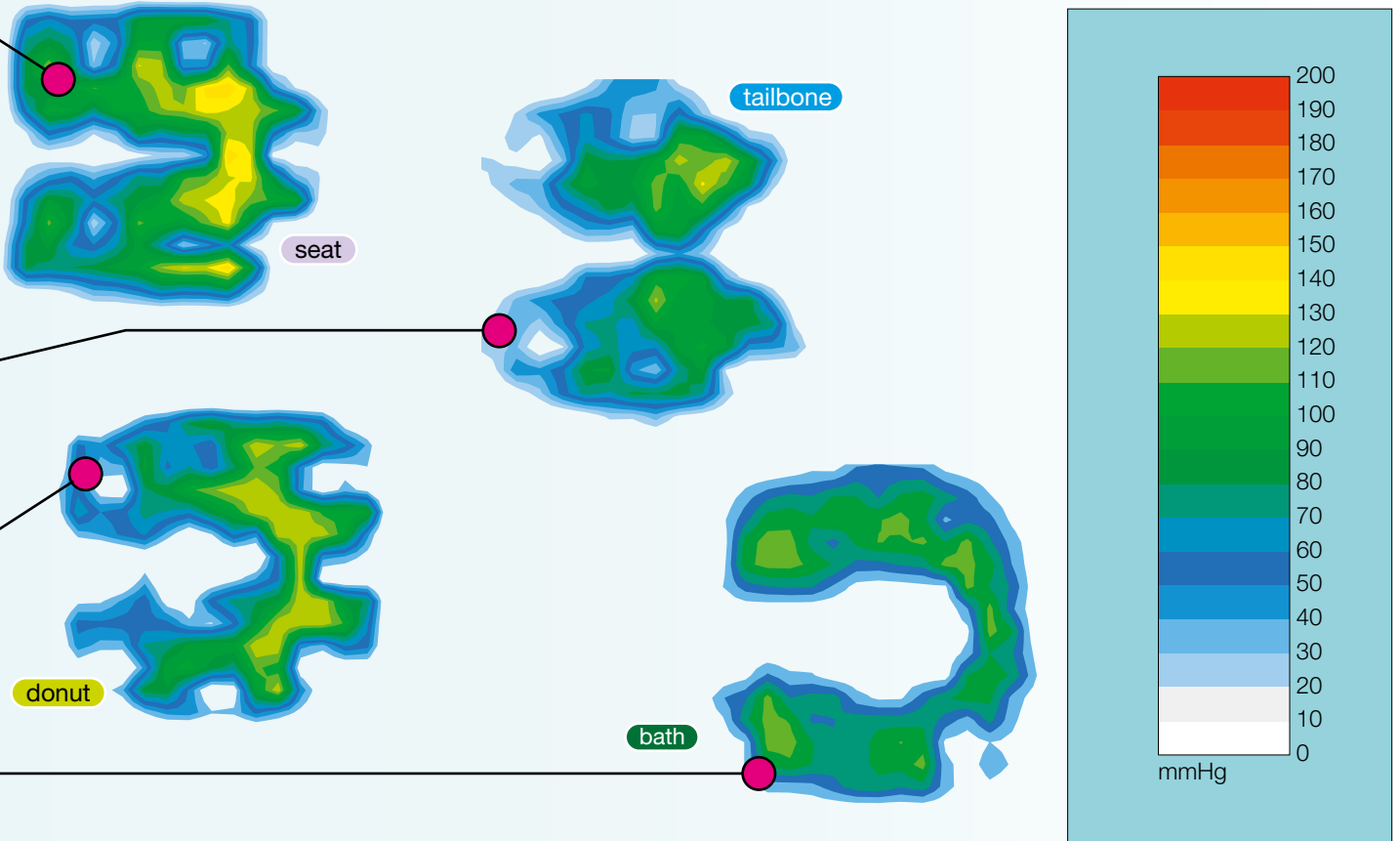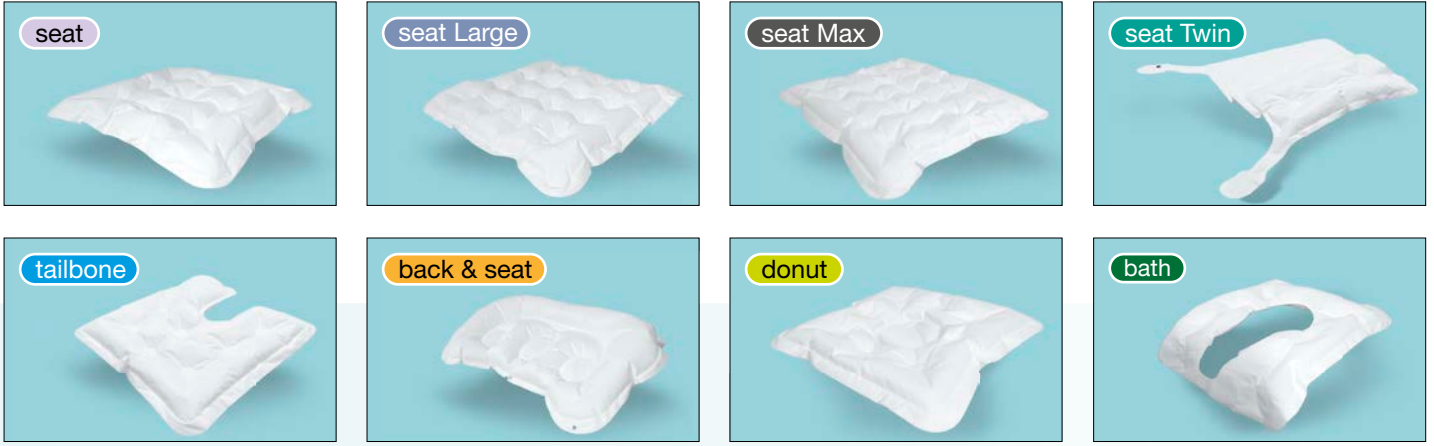

Mattress with differentiated pressure relief






Multi chamber technology

all on[®]

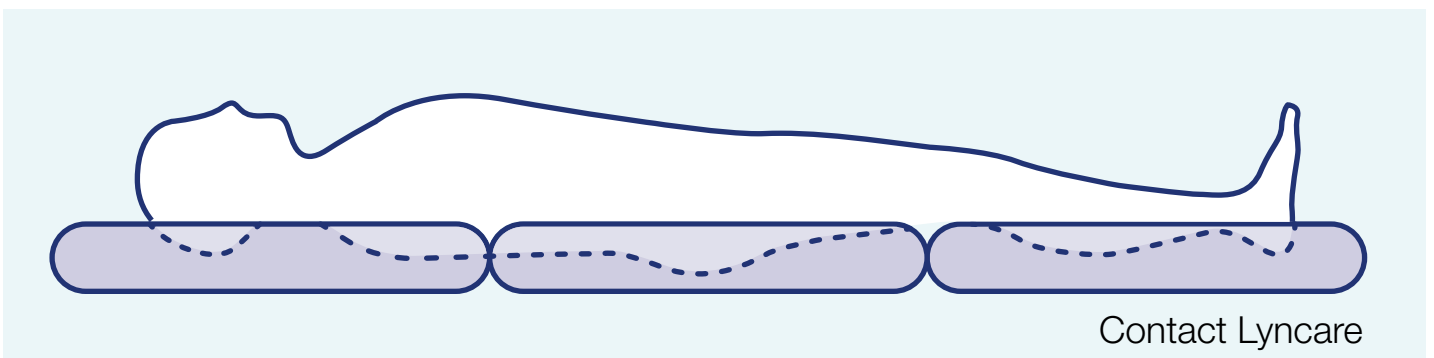

Cushions for pressure relief of body and arms

Name	Size in cm.	Item no.	About the device	Application
All up - Multi Small	50 x 25 x 8	70003	Positioning cushion with 2 chambers. Multifunctional. Can be used for both supine and sitting positions. (between legs, underneath the back, providing excellent support under arms and between the knees.).	
All up - Multi Medium	50 x 37 x 8	70004	Positioning cushion with 3 chambers. Multifunctional. Especially useful for a supine position (between legs, back), and possibly as support under the feet.	
All up - Multi Large	50 x 50 x 8	70005	Positioning cushion with 4 chambers. Multifunctional. Can be folded/rolled up. For use under the leg or arm, comes with an extra Velcro strap.	
All up - Multi Max	50 x 50 x 20	70012	Positioning cushion, particularly useful for hip fracture patients and when the lower limb needs to be elevated above heart level.	
All up - Universal	50 x 50 x 8	70010	Positioning cushion, typically for the back, sacrum and tailbone in adjustable bed and the lateral positions between the legs. Can be folded if a higher cushion is desired.	
All up - Universal Long	55 x 70 x 8	70015	Positioning cushion, typically for the lateral positions between the legs, longer model which protects the knees and malleolus. Can be folded in half if a higher cushion is desired.	
All up - Universal Slim	45 x 75 x 8	70016	Positioning cushion, typically for lateral positions between the legs, longer and narrower model which ensures a distance between the knees in the lateral position. Can be folded if a higher cushion is desired.	
All up - Head	28 x 30 x 5	70011	Cushion, particularly useful for the preventing and treating of pressure sores on the ears and neck, and is also useful for coma patients.	
All up - Wing	50 x 85 x 8	70006	Relief cushion for lower limbs, used for bed-bound patients. The cushion is easily fixed to the mattress.	
All up - Wedge 30 deg.	70 x 40 x 15	70024	Positioning cushion for the lateral position. Ensures a stable position at 30°, and the wedge shape is fixed by a flap that is pushed underneath the patient.	
All up - Lap	80 x 60 x 10	70025	Positioning cushion, crescent. Advantageous for wheelchair users for supporting the arms.	
All up - Lap Max	80 x 60 x 20	70026	Positioning cushion, crescent. With extra height. Advantageous for wheelchair users for supporting the arms.	
All up - Arm Rest	55 x 25 x 8	70027	Positioning cushion for supporting the arms.	
All up - Elbow	30 x 20 x 10	70028	Positioning cushion for supporting the elbow.	

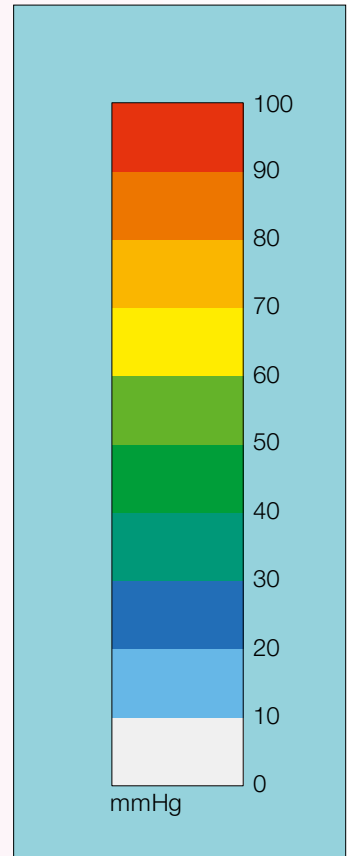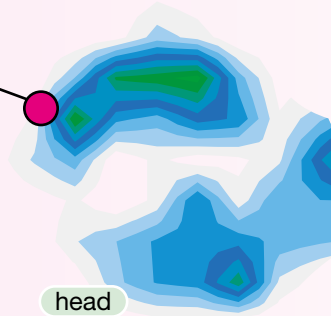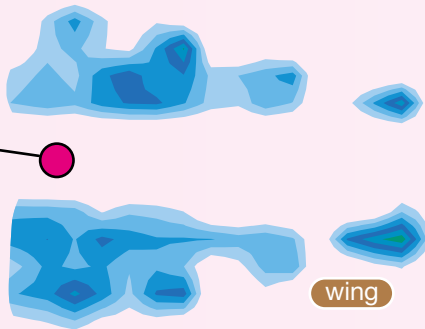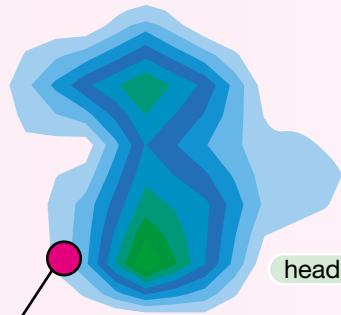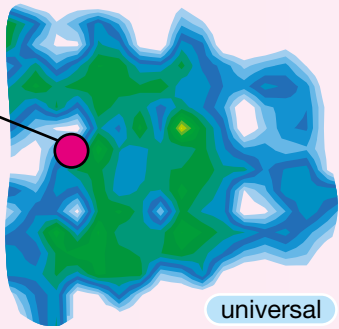

Levabo Emergency Package

pressure ulcer 
dcubi kit[®]

Our products really make a difference.

The innovative positioning products from Levabo[®] are developed in collaboration with leading wound care nurses and are produced in Denmark. The products have a series of advantages compared to conventional positioning products.

With the use of Levabo's positioning cushions you always have the required device at hand for effective prevention of pressure sores. The products are single-patient, inflatable positioning cushions, which can be disposed of after use. The products take up very little space and are easily transportable – as is Levabo's emergency package, developed for this very purpose.

These advantages are especially useful for the domiciliary care sector. As a domiciliary care nurse you can initiate effective prevention of pressure sores at an early stage with Levabo's products, much easier and cheaper than with conventional positioning products. Immediate action can be taken to position the patient optimally, and thereby avoid unnecessary discomfort for the patient, in addition to saving the health care system potentially expensive future treatment.

The products are multifunctional and can be used in several different supportive ways for patients. They are single-patient use, thereby ensuring optimum levels of hygiene.



Name	Size in cm.	Item no.
dcubi kit bag	40 x 31 x 15	76007

Levabo Guide


For your iPhone, iPad & Android



Get the free Levabo app on:

<http://www.levabo.dk/apps.html>

Or scan the QR code:

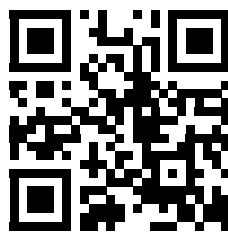

Levabo 
academy

visit www.levabo.com



Lyncare, CB, 28.08.2017

Distributed in Ireland by



Lyncare Ltd, Unit 3 Shannonside Business Park,
Birdhill, Co. Tipperary, V94 V8W3
Tel: 061 379090, Fax: 061 379080,
Email: info@lyncare.ie, Web: www.lyncare.ie

LEVABO[®]

medical cushions



Levabo ApS
Højlundsparken 49
8355 Solbjerg
Tel.: +45 20 70 09 30
info@levabo.dk
www.levabo.com